


Product	PLUS shower seat 450, manually height adjustable	
Item number	<p>R7430182112 VVS: 777535912 RSK: 7973585 EAN: 5708590312002 HMI: HMI 71800 NBS: N13:448 FDA: INN/D135865</p>	
Colour	Anthracite grey	
Description	<p>Shower seat 450, with back- and armrests. Manually height adjustable by 200 mm. Armrests are independently height adjustable. The seat and the armrests are foldable. Pale grey cover plate for vertical wall track.</p> <p>Materials: Aluminium and the shower seat consists of a plywood core covered in comfortable polyurethane foam. The non-slip gripping surface of the armrests are made of soft SEBS incorporating a glass-reinforced nylon core.</p> <p>Excl. fixing kit.</p>	
Adjustability	Height adjustable	
Adjustment height	200 mm	
Maximum load	Maximum load: Seat: 300 kg, armrest: 80 kg.	
Material	<p>The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with an in-mould coated surface that is protected using a clear PU lacquer. The backrest is made of polyurethane (PU) foam protected using a clear PU lacquer. The arms are primarily constructed of 10-micrometre anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.</p>	
Ergonomy	<p>Ergonomically designed seat - contact surface are 450 x 437 mm. Soft and slip-resistant material that has a warm feel. The oval shape of the armrests gives an easy and comfortable surface to lean on and push against.</p>	



Product	PLUS shower seat 450, manually height adjustable
Features	The seat, backrest and armrests can be adjusted simultaneously thanks to a safe and easy click system. The armrests can also be adjusted independently of the seat and backrest. The seat and armrests can fold up when not in use or to maximise floor space.
Cleaning	<p>Pressalit Care recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in clear water.</p> <p>This product complies with the requirements of Council Directive 93/42/EEC concerning medical devices, as last amended by Council Directive 2007/47/EC of 5 September 2007. It also complies with all relevant requirements in the Council Directive 2006/42/EC (of 17 May 2006) on machinery.</p>

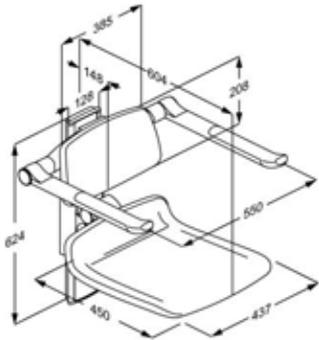

Schrama Handels GmbH
Schwendberg 465
A-6283 Hippach, Austria
+43 5285-63067
www.hygieneprofi.biz