


Product	PLUS shower seat 450, electrically height adjustable	
Item number	<p>R7470112  VVS: 777537012  RSK: 7973686  EAN: 5708590313443  HMI: HMI 71816  NBS: N13:448</p>	
Colour	Anthracite grey	
Description	<p>Powered shower seat 450, with back- and armrests. Electrically height adjustable by 500 mm, the armrests are independently height adjustable. The seat and the armrests are foldable.</p> <p>Materials: Aluminium and the shower seat consists of a plywood core covered in comfortable polyurethane foam. The non-slip gripping surface of the armrests are made of soft SEBS.</p> <p>Incl. fixing kit for brick wall.</p>	
Adjustability	Height adjustable	
Adjustment height	500 mm	
Maximum load	Maximum load: 300 kg, armrests 80 kg. During height adjustment, maximum load is 150 kg.	
Speed without load	Speed up: 18 mm/s Speed down: 25 mm/s	
Noise level	45 dB	
Material	<p>The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with an in-mould coated surface that is protected using a clear PU lacquer. The backrest is made of polyurethane (PU) foam protected using a clear PU lacquer. The arms are primarily constructed of 10-micrometre anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.</p>	



Product	PLUS shower seat 450, electrically height adjustable
Ergonomy	Ergonomically designed seat - contact surface are 450 x 437 mm. Soft and slip-resistant material that has a warm feel. The oval shape of the armrests gives an easy and comfortable surface to lean on and push against.
Features	Powered raising or lowering at the push of a button - the seat can easily be adjusted upwards or downwards to a good working height for the carer while the user sits on it. The seat, backrest and armrests can be adjusted simultaneously thanks to a safe and easy click system. The armrests can also be adjusted independently of the seat and backrest. The seat and armrests can fold up when not in use or to maximise floor space. Wired hand control comes with brackets for fitting to the wall.
Mains connection (V)	230
Actuator (motor)	LA27CS (Linak) 6000N, IPX6
Control box	TR6 (Linak) IPX0
Hand control	HB7 (Linak) IP66
Power usage	Power consumption by maximum load: 150 W Power consumption by standby: 2 W
Cleaning	<p>Pressalit Care recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in clear water.</p> <p>This product complies with the requirements of Council Directive 93/42/EEC concerning medical devices, as last amended by Council Directive 2007/47/EC of 5 September 2007. It also complies with all relevant requirements in the Council Directive 2006/42/EC (of 17 May 2006) on machinery.</p>

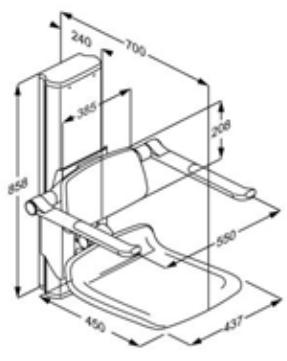

Schrama Handels GmbH  
 Schwendberg 465  
 A-6283 Hippach, Austria  
 +43 5285-63067  
 www.hygieneprofi.biz