


| | | |
|--------------------|---|---|
| Product | Nursing bench 3000, electrically height adjustable | |
| Item number | <p>R8513021 VVS: 777540234 EAN: 5708590265575 HMI: HMI 42881 NBS: N13:40,375</p> |  |
| Colour | Lead grey | |
| Description | <p>Nursing bench with electric motor, 1300 mm, including water collection tray with integrated outlet Height adjustable 700 mm</p> <p>Incl. fixing kit for brick wall</p> | |
| Adjustability | Height adjustable | |
| Adjustment height | 700 mm | |
| Maximum load | Maximum load: 200 kg. | |
| Speed without load | <p>Speed up: 15 mm/s Speed down: 20 mm/s</p> | |
| Noise level | 48 dB | |
| Material | Aluminium, polyurethane foam, plastics, powder-laquered steel | |
| Ergonomy | <p>Pressalit Care offers a range of high quality, fixed or power adjusted nursing benches designed to provide a stable platform suitable for changing adults and children in a variety of settings: schools, day care centres, institutions, hospitals and private homes. The nursing bench 3000 series provides for maximum safety and comfort. Besides an award winning design, the series includes comfortable and durable surfaces, including polyurethane foam and aluminum.</p> <p>The flexibility of the nursing bench provides for optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p> | |

| | |
|---|--|
| Product | Nursing bench 3000, electrically height adjustable |
| <p>Features</p> <p>Mains connection (V)</p> <p>Actuator (motor)</p> <p>Control box</p> <p>Hand control</p> <p>Power usage</p> | <p>A water collection tray with integrated water outlet prevents water from splashing on carers or the floor. The flexible hose attached to the centre mounted outlet of the water collection tray, can be connected to either an outlet in the floor or on the wall. Adjustable neck support which can be moved to either end of the 1800 mm nursing benches. The slightly curved form makes the bench extremely comfortable and safe. The integrated gas cylinders enable the bench to be folded out of the way when not in use to maximise floor space.</p> <p>230</p> <p>LA31 Careline, (Linak), 6000 N, IPX6</p> <p>CB9 Careline, (Linak), IPX6</p> <p>HB7 (Linak), IP66</p> <p>Power consumption by maximum load: 70 W Power consumption by standby: 2 W</p> |
| Cleaning | <p>Pressalit Care recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in clear water.</p> <p>This product complies with the requirements of Council Directive 93/42/EEC concerning medical devices, as last amended by Council Directive 2007/47/EC of 5 September 2007. It also complies with all relevant requirements in the Council Directive 2006/42/EC (of 17 May 2006) on machinery.</p> |

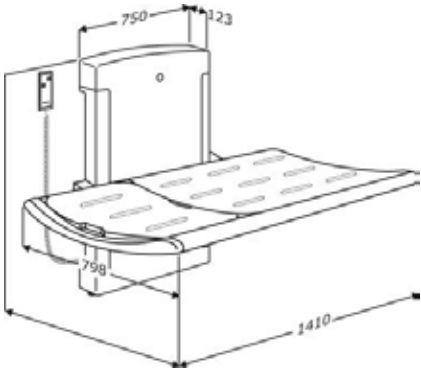